*[®] nosiboo[®]pro

electric nasal aspirator

reddot award 2014
winner

Instructions for Use

EN

AR

BG

CS

DA

DE

EL

ES

ET

FI

FR

HR

HU

IT

LT

LV

NL

NO

PL

PT

RO

SK

SL

SV

TR

VI

ZH

Table of contents

EN	Instructions for Use	LT	Naudotojo vadovas
AR	دليل المستخدم	LV	Lietošanas pamācība
BG	Ръководство за употреба	NL	Gebruiksaanwijzing
CS	Návod k použití	NO	Bruksanvisning
DA	Brugervejledning	PL	Instrukcja obsługi
DE	Bedienungsanleitung	PT	Manual de usuário
EL	Εγχειρίδιο	RO	Instrucțiuni de utilizare
ES	Manual del Usuario	SK	Návod na použitie
ET	Kasutusjuhend	SL	Navodila za uporabo
FI	Käyttöoppaan	SV	Bruksanvisning
FR	Manuel d'utilisation	TR	Kullanım Kilavuzu
HR	Upute za uporabu	VI	HƯỚNG DẪN SỬ DỤNG
HU	Használati útmutató	ZH	使用说明书
IT	Manuale d'istruzioni		

On/off
control button

Pro device

Power cord

Air mesh (solid objects will
not enter the housing)

Connecting unit

Tube


Colibri head

Colibri head part without
connector

Collared tube

Colibri head part with
connector

Nose tip

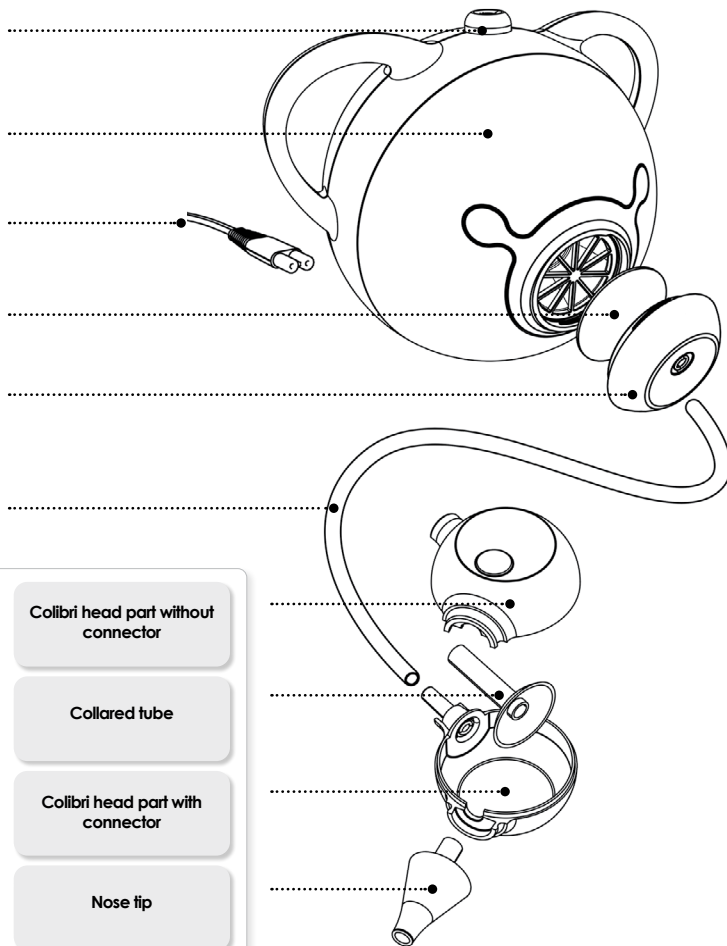Instructions for Use

Foreword

Dear Customer,
thank you for purchasing our KSCI, PTPA, Red Dot Design and German Design award-winning Nosiboo Pro nasal aspirator.

Please, follow the instructions provided in this Instructions for Use to obtain the very best performance from your nasal aspirator.

Table of contents

1. Package contents
2. Safety instructions
3. Intended use of the device
4. Before use
5. Instructions for use
6. Cleaning
7. Storage and maintenance
8. Troubleshooting
9. Important information about electromagnetic interference
10. Environmental impact
11. Specifications
12. Accessories, spare parts
13. Warranty

Read all instructions before use.

CE 2409

Nosiboo Pro electric nasal aspirator is a medical device. You may exercise your right of withdrawal only if the device has not been used, is in undamaged condition, and the hygienic packaging of the device, the Colibri head, and the nose tip remains unopened. In all other cases, the purchase price cannot be refunded.

1. Package contents

Nosiboo Pro

Tube

Colibri head +
Newborn's nose tip

Standard nose tip

Power cord

Cleaning brush

ATTRACT Kft. is an ISO 13485 certified medical device manufacturer. This product is an MDD Legacy Device.

EUDAMED SRN: HU-MF-000004252

2. Safety instructions

Symbols

Class IIa medical device

Read the Instructions for Use

GTIN – Global Trade Item Number (UDI-DI)

Lot number

Serial number

Model number (catalogue number)

Unique Device Identifier

Protection against solid particles of diameter greater than 12.5 mm; protected from condensation.

Double insulation

Name and address of the manufacturer

Keep out of reach of children

Keep away from sunlight

Fragile – Handle with care

Keep dry

Collect separately the electrical and electronic equipment waste!

Caution!

General use

Use only as directed. Misuse can cause injury or death.

The device should only be operated by adults.

Do not allow children to play with the device. It is NOT A TOY. (The power cord can easily wind around a child's neck.)

Do not leave the device unattended near children.

Do not use the device in the bathroom or in any other humid or wet environment.

Do not handle the power cord with wet hands.

Do not cover the vents on the back of the device while in use. Hot air exits the device through the vents.

Do not use the nasal aspirator for nasal washing!

The device should be on a flat, stable surface during use.

The connection point of the Colibri head and the tube should point downwards while in use.

Always switch off the device after use by turning the control button. Do not switch off the device by unplugging the power cord!

Always thoroughly clean the Colibri head after use (see the "Cleaning" section).

Read and follow the instructions in the "Important information about electromagnetic interference" section.

Read and follow the instructions in the "Environmental impact" section if you would like to dispose of the device or any part of it.

The device is intended to be used only for removing nasal mucus from a child's nose.

Please consult your healthcare provider about when and how often you should clear your baby's nose, or if redness occurs.

Keep out of reach of children.

Make sure that all parts are properly attached before using the device.

Make sure that the housing, the tube, and the Colibri head are not damaged or cracked.

Stop using the device if mucus or saline solution has been drawn into the tube. Disconnect it from the electrical outlet, empty the Colibri head, clean and dry the Colibri head and the tube, check the air mesh and replace it if it is wet or clogged.

Double-check that you have inserted the collared tube into the Colibri head after cleaning.

Do not insert the nose tip into the baby's mouth or ears.

Do not force or insert the nose tip too far into the nose.

Do not press the nose tip against the inner side of the nostril.

Do not use the device if nasal sores, inflammation, or nosebleed are present. If pain or bleeding occurs, stop use immediately and consult your doctor.

Do not use the device if you see that the tube or the Colibri head are damaged. Replace the damaged part. (A list of spare parts can be found in the "Accessories, spare parts" section.)

Do not use the Colibri head and the device without the collared tube.

Do not use more than one teaspoon of saline solution. Excess solution can enter the connecting tube and can cause damage to the device.

The device cannot be used in oxygen-rich environments and in clinical circumstances! The medical device is for home use only.

Do not attach the Colibri head to other devices.

Do not unplug the power cord from the electrical outlet until the device is completely shut down.

Cleaning, Maintenance, and Storage

Protect the power cord against damage.

Using a damaged power cord can cause short circuit (short-circuit fire).

Keep the power cord away from heat (e.g., furnace, heater, or convector) because it can cause damage to the insulation.

Turn the device off and grasp the plug to unplug the power cord from the electrical outlet when assembling, disassembling, and cleaning the device.

Do not break the power cord, and do not place heavy items on the cord, to avoid damage. Grasp the plug to unplug the power cord from the electrical socket. A damaged cord may cause electric shock.

Do not remove the mounting screws from the housing. Do not reconstruct the device for other purposes.

A broken, damaged, worn, frayed or short-circuited power cord can cause an electrical fire, so do not use it!

Do not leave the parts of the dismounted Colibri head within reach of children. Children may swallow the small parts.

Do not wash the device housing or the power cord in water or any other liquid. Do not handle the power cord with wet hands. This may cause electric shock or malfunction.

Make sure that all parts are dried thoroughly before the next use.

Double-check that you have inserted the collared tube into the Colibri head after cleaning.

When not in use, always unplug the device from the outlet by grasping the plug.

To prevent damage, the connecting unit should not be exposed to heat.

It is not allowed to use a dishwasher, boiling water, microwave oven, vaporizer, or baby bottle sterilizer for cleaning the spare parts, because they might deform by heat!

Store in a cool, dry place. Avoid direct sunlight or heat.

Power Source

Protect the power cord against damage.

Turn the device off and grasp the plug to unplug the power cord from the electrical outlet when assembling, disassembling, and cleaning the device.

In case of unexpected failure or malfunction, immediately switch off the device and unplug the power cord from the electrical outlet.

Be sure to use the device only with a 230 V 50 Hz AC.

Do not use the device if there is a breakage or crack in the housing or power cord.

Do not break the power cord, and do not place heavy items on the cord, to avoid damage. Grasp the plug to unplug the power cord from the electrical socket. A damaged cord may cause electric shock.

When not in use, always unplug the device from the outlet by grasping the plug.

Do not unplug the power cord from the electrical outlet until the device is completely shut down.

3. Intended use of the device

The Nosiboo Pro electric nasal aspirator is a medical device intended for use by laypersons. It is designed for use at home to remove nasal mucus from children primarily under the age of 6. Its purpose is to help manage symptoms of cold and nasal congestion, and to prevent sinus infections in children.

The device is not intended for use by children and persons with reduced physical, sensory or mental capabilities.

Intended purpose:

- * to ensure clearer nasal breathing for infants,
- * to relieve symptoms of common childhood colds.

Do not use:

- * if the child has haemophilia,
- * if the child has frequent nosebleeds,
- * if the child is taking blood-thinning medicines,
- * in the event of congenital malformations (e.g. cleft palate).


The device is not a nasal irrigator! Do not use large amounts of liquid during nasal aspiration! The Colibri head can hold up to a teaspoon (5 ml) of thin liquid.

4. Before use

1. Disassemble the Colibri head and clean and dry its parts.

Parts of the Colibri head:

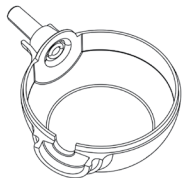
Colibri head part with connector

Colibri head part without connector


Nose tip

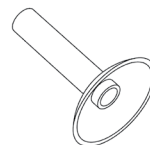
Collared tube


2. Assemble the Colibri head.

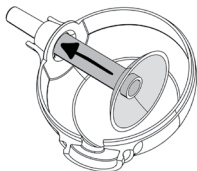
a) Insert the collared tube into the side of the Colibri head part which is supplied with a connector.

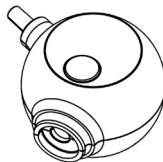b) Place the Colibri head part that has no connector in the correct position.

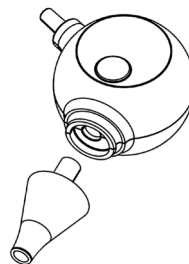
c) Push the nose tip into place.

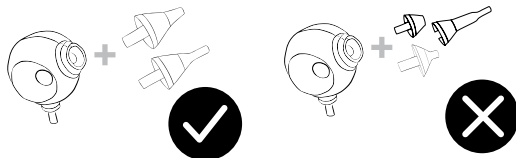
The use of non-original spare parts poses technical and medical risk for which the manufacturer accepts no responsibility! The use of non-original spare parts may cause malfunction or damage to the device and will result in loss of warranty. Only use official Nosiboo parts!

3. Attach the Colibri head to the device.

- Attach one end of the connecting tube to the Colibri head.
- Attach the free end of the connecting tube to the device housing.

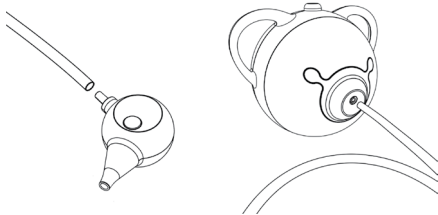4. Attach one end of the power cord to the device.

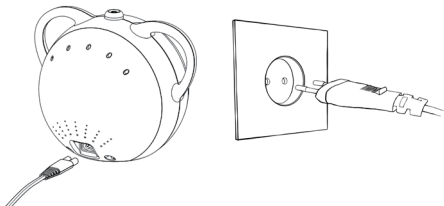Plug the other end of the power cord into a 230 V 50 Hz AC outlet.

5. The device is ready for use.

The device should be connected to or disconnected from the power network only if the control button is in OFF position (the device is switched off).

Be sure to use the device only with a 230 V 50 Hz AC. Failure to do so may cause electric shock or fire.

Clean the parts of the Colibri head before every use.

Make sure that all parts are clean and dried thoroughly before use.


Do not use the Colibri head and the device without the collared tube.

5. Instructions for use

Do not aspirate water with the device!

Do not use the nasal aspirator for nasal washing!

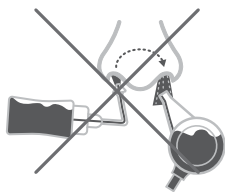Aspirating water with the electric device will result in loss of warranty.

1. The device is ready for use.

2. Switch on the device by turning the control button.


The device should be connected to or disconnected from the power network only if the control button is in OFF position (the device is switched off).

NOTE: Any time you are not actively using the device, switch it off. After the motor cycles off completely, you can restart it.

3. Select the correct suction power.

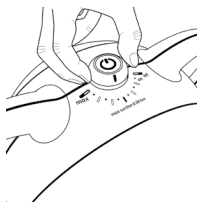low setting for newborns

medium setting for babies

high setting for small children and for preschoolers

TIP: If the mucus is thin, you can use the device at a lower setting. If it is thick, you can use the device at a higher setting than normal.

4. Hold your child in a stable position.

Sit your child onto your lap, hold them in your arms or sit them into an upright position supporting their back with a pillow. You may also sit them in an infant bouncer or infant seat.

5. Insert the clean, dry nose tip into the child's nostril.

In case of newborns, it is sufficient to hold the nose tip at the opening of the nostril. The connection point of the connecting tube should point downwards while in use.

NOTE: In case of swollen nasal mucosa or blocked nasal openings, use decongestant solutions or nasal sprays 30 minutes before application.

6. Move the nose tip around gently to make sure that you get all areas of the nostril.

NOTE: Do not pinch or block the other nostril while using the nasal aspirator.

Newborn's nose tip:
0-6 months (baby's age)

Standard nose tip:
6 months and older


The suggested duration of contact aspiration is 10 seconds per nostril.

Dry, very thick mucus can be dissolved by a mild saline solution. After using the solution, wait a few minutes, wipe excess discharge with a tissue, then reuse the nasal aspirator as directed.

To avoid damage to the device, do not use more than a teaspoon of solution at a time. The Colibri head can hold up to a teaspoon (5 ml) of thin liquid.

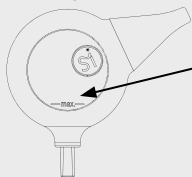When the Colibri head fills to the marked level (max.), empty it to prevent the secretion from entering the tube.

This device is not a nasal irrigator.

The device should only be operated by adults.

Do not allow children to play with the device. It is NOT A TOY.

Do not cover the vents on the back of the device while in use. Hot air exits the device through the vents.

Do not use the device in the bathroom or in any other humid or wet environment.

Do not handle the power cord with wet hands.

Do not break the power cord during use and storage.

7. Switch off the device after use with the control button.

The device should be connected to or disconnected from the power network only if the control button is in OFF position (the device is switched off).

NOTE: Always switch off the device after use.

TIP: The Nosiboo Pro Nasal Aspirator can be used several times a day. It is recommended that you clear the airways before bedtime and feeding to ensure the child breathes peacefully.

NOTE: Nose drops and other preparations should be used as directed by your healthcare professional.

Make sure that all parts are clean and dried thoroughly before use.

Always switch off the device after use by turning the control button. Do not switch off the device by unplugging the power cord!

The device should be on a flat, stable surface during use.


Do not unplug the power cord from the electrical outlet until the device is completely shut down.

Do not use the Colibri head and the device without the collared tube.

Do not cover the device, and do not use it on a soft surface which could block the ventilation holes (couch, pillow, bedclothes etc.)

6. Cleaning

1. Remove the connecting tube and the nose tip, then take apart the Colibri head.

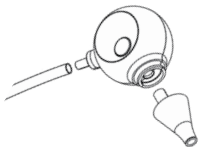
2. Wash the collared tube, the nose tip, and the Colibri housing in lukewarm water and dishwashing detergent and let them dry.


Use the brush for cleaning the nose tip.

3. Clean the housing with a dry, soft cloth; do not use any liquids, detergents, or chemicals.

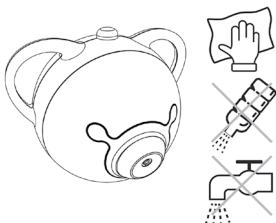4. Wipe the connecting tube with a cloth if necessary.

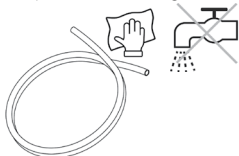Turn the device off and grasp the plug to unplug the power cord from the electrical outlet when assembling, disassembling, and cleaning the device.

Do not wash the device housing or the power cord in water or any other liquid. Do not handle the power cord with wet hands. This may cause electric shock or malfunction.

Do not wash the connecting tube.

Thoroughly clean the Colibri head after use.

Make sure that all parts are clean and dried thoroughly before use.

If liquid gets into the connecting tube due to improper use, make sure that it is completely dry before using it again.

Do not remove the mounting screws from the housing. Do not reconstruct the device for other purposes.

To avoid damage, make sure the device is not exposed to heat.

It is not allowed to use a dishwasher, boiling water, microwave oven, vaporizer, or baby bottle sterilizer for cleaning the spare parts, because they might deform by heat!

Do not use any liquids, detergents or chemicals for cleaning the device housing and the connecting tube.


7. Storage and maintenance

Storage

1. Unplug the power cord from the wall socket only after the device has completely shut off.

NOTE: After the device is switched off, the power cord should be unplugged from the electrical outlet and from the device by grasping the plug.

2. Disconnect the connecting tube from the device and the Colibri head.

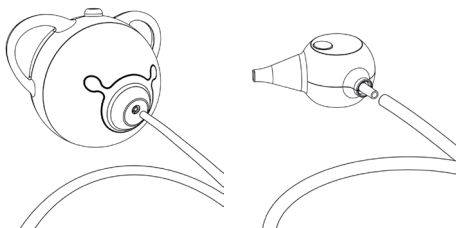3. Keep the device in an enclosed space where it is protected from dust and damage (e.g., in a cabinet).

NOTE: Do not store near a heat source.

Maintenance

There is a special air mesh behind the connecting unit which prevents solid objects and dust from entering the housing. In the event of an unexpected drop in performance, it is recommended that you check the air mesh. If the air mesh is dirty or clogged, replace it with a spare air mesh. (A list of spare parts can be found in the "Accessories, spare parts" section.)

Protect the power cord against damage.

Keep the power cord away from heat.

Do not leave the parts of the dismantled Colibri head within reach of children. Children may swallow the small parts.

Do not break the power cord, and do not place heavy items on the cord, to avoid damage. Grasp the plug to unplug the power cord from the electrical socket. A damaged cord may cause electric shock.

Do not wind the power cord around the device.

A damaged power cord can cause electric fire!

Make sure that all parts are clean and dried thoroughly before use.

Store the device in a clean, dry place. Avoid direct sunlight or heat.

8. Troubleshooting

Malfunction	Reason	Recommendation
The device cannot be switched on.	The device is not connected to the power network.	Connect the power cord to the device, then to the electric socket.
The suction intensity is decreased.	The nose tip is blocked or the Colibri head was incorrectly assembled.	Take the Colibri head apart and clean it. After drying the parts of the Colibri head, assemble them according to the Instructions for Use.
The suction intensity is decreased.	The air mesh is clogged.	Replace the air mesh according to the Instructions for Use's "Storage and maintenance" section.
Liquid or mucus got into the connecting tube.	The Colibri head was not held in the correct position.	Switch off the device and clean and dry the Colibri head and the connecting tube. Check the air mesh and replace it if necessary.
Liquid or mucus got into the connecting tube.	The connection point of the tube is not pointed downwards.	Switch off the device and clean and dry the Colibri head and the connecting tube. Check the air mesh and replace it if necessary.
Liquid or mucus got into the connecting tube.	Too much saline solution was used.	Reduce the volume of the saline solution to less than a teaspoon (5 ml). Switch off the device and clean it, then dry the Colibri head and the connecting tube. Check the air mesh and replace it if necessary.

9. Important information about electromagnetic interference

Home wireless network equipment, mobile phones, and wireless phones may affect the operation of certain medical equipment, as electromagnetic interference can occur between devices. Electromagnetic interference may cause harmful malfunction in the operation of certain medical devices, potentially causing emergency situations. It is recommended to keep at least 5 m distance between the devices!

Due to the abovementioned reasons, medical devices must not interfere with any other equipment. Do not use near high-performance radiation-emitting devices.

The Nosiboo Pro Nasal Aspirator is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled.

The Nosiboo Pro nasal aspirator is certified as compliant with the requirements of IEC 60601-1-2:2014; EN 6061-1-2:2018 (EMC test); EN 6061-1-11:2015; EN 6061-1:2017 standards.

10. Environmental impact

The device should not be disposed of in a household waste bin but should be disposed of as electrical waste in accordance with local regulations. Take it to the designated collection point.

The Colibri head does not present a risk to the environment and can be disposed of as household waste.

11. Specifications

Max. motor performance: 500 W
Adjustable motor performance: 150-390W
Performance: 420 VA
Differential pressure: 0.04-0.08 bar
Max. air flow: 30-42 l/min
Power supply: 230 V 50 Hz
IP class: IP21
Protection class: II
Operation type: no continuous operation
Operation length: on 15 min/appr. off 50 min
Duration of continuous use: max 15 min.
Safety switch release time: after 15 minutes
Tripping temperature of the thermal switch: 95 °C
Average housing temperature: 30-35 °C after 1 min continuous use
Noise level: 60-70 dB
Expected life of the housing (electrical parts): 3 years
Expected life of the Colibri head: 3 months
Diameter of the nasal insertion end of the standard nose tip:
max. 7.0 mm
Length of the nasal insertion end of the standard nose tip:
max. 8.0 mm
Diameter of the nasal insertion end of the newborn's nose tip:
max. 4.9 mm
Length of the nasal insertion end of the newborn's nose tip:
max. 12.6 mm
Weight: 1090 g
Dimensions (housing, W x H x D): 198 x 162 x 172 mm
Dimensions (complete package): 180 x 180 x 225 mm
Operating conditions:
1. Temperature range: +5 °C to +40 °C
2. Relative humidity: max. 75%
3. Atmospheric pressure / barometric pressure range: 525-795 mmHg / 70-106 kPa
Transportation conditions:
1. Temperature range: -20 °C to +65 °C
2. Relative humidity: 10-75%
3. Atmospheric pressure / barometric pressure range: 150-795 mmHg / 20-106 kPa
Storage conditions:
1. Temperature range: 0 °C to +40 °C
2. Relative humidity: 10-75%
3. Atmospheric pressure / barometric pressure range: 525-795 mmHg / 70-106 kPa
MD Category: II.a
MDD/MDR: MDD Legacy Device
EUDAMED SRN: HU-MF-00004252

We reserve the right to make technical changes.

12. Accessories, spare parts

In the European Union, parts can be ordered from the specialist service station:

ATTRACT Kft. H-7622 Pécs, Siklósi út 1/1. - Hungary
Tel.: +36 72 551 642 (Mo-Fr: 8 am – 4 pm)
E-mail: customerservice@nosiboo.com

In countries outside the European Union, spare parts are available at the distributors of Nosiboo products.

Power cord EU
300026

Power cord UK
300054

Connecting unit
300007 * blue
300008 * pink
300009 * green
300076 * grey
300077 * sage
300122 * beige

Nosiboo Pro Accessory Set:

Accessory pack consisting of Nosiboo Pro/Pro2 spare parts. All the parts included in the package can be used for both types of the device.

Available at the manufacturer's webshop (<https://shop.nosiboo.eu/>), at Amazon, or the distributors of Nosiboo products.

- * Colibri head with newborn's nose tip
- * standard nose tip
- * tube
- * air mesh pack
- * cleaning brush

Available colours:

- * ACS-01-01 blue
- * ACS-01-02 pink
- * ACS-01-04 grey
- * ACS-01-05 beige
- * ACS-01-13 sage

Manufacturer contact details

ATTRACT Kft.

H-7622 Pécs, Siklósi út 1/1., Hungary

Tel.: +36 72 551 642 (Mo-Fr: 8 a.m. – 4 p.m.)

E-mail: contact@nosiboo.com

13. Warranty

The warranty period provided by the manufacturer for the device from the date of purchase is: 24 months.

The warranty does not apply to products used improperly, and the manufacturer cannot be held liable for any damage resulting from improper use. The manufacturer is not responsible for defects in products originating from unofficial import sources.

If the warranty conditions are met, the distributor is responsible for handling the warranty claim. The distributor must repair or replace the device without undue delay. The distributor is the issuer of the invoice confirming your purchase. Please keep the invoice for 24 months for warranty purposes.

For products placed on the market in the European Union:

Please keep the warranty card and the purchase invoice for 24 months for warranty purposes. The warranty card included in the device packaging is valid only together with the invoice. For more information about the warranty, please visit the link below:
<https://nosiboo.eu/downloads>

What is not covered by the warranty?

- * Improper use.
- * Products that have not been cleaned or maintained properly.
- * Mechanical damage or disassembly of the device housing. Disassembly of the housing for any purpose other than battery replacement will void the warranty. Damage visible on the device housing or bottom label may indicate disassembly.
- * Repairs or repair attempts not performed by an authorized service centre.

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